

FIRST AID KIT&MEDICAL KIT - INSPECTION/CHECK

EFFECTIVITY: ACFT MODEL(S) EMB-135

1. General

- A. This section gives the procedures for the inspection/check of the first-aid kit and medical kit.
- B. The first-aid kit and medical kit are installed on the LH wardrobe, RH wardrobe, outboard galley, forward galley or aft galley, as applicable.
- C. On the aircraft equipped with 16 seat, the first-aid kit and medical kit are installed on the LH Forward Stowage.
- D. The procedures in this section are given in the sequence below. The tasks identified with (♦) are part of the Scheduled Maintenance Requirements Document (SMRD).

TASK NUMBER	DESCRIPTION	EFFECTIVITY
25-60-01-200-801-A ♦	FIRST-AID KIT AND MEDICAL KIT - INSPECTION	ACFT MODEL(S) EMB-135

TASK 25-60-01-200-801-A

EFFECTIVITY: ACFT MODEL(S) EMB-135

2. FIRST-AID KIT AND MEDICAL KIT - INSPECTION

A. General

- (1) The function of this inspection is to make sure that the condition of the first-aid kit and medical kit is good and that their permitted life is not expired.

B. Zones and Accesses

Not Applicable

C. Tools and Equipment

Not Applicable

D. Auxiliary Items

Not Applicable

E. Consumable Materials

Not Applicable

F. Expandable Parts

Not Applicable

G. Persons Recommended

QTY	FUNCTION	PLACE
1	Does the task	On the wardrobe shelf

H. Inspect (General Visual) First-Aid Kit and Medical Kit for Condition and Validation

SUBTASK 212-005-B

- (1) Remove the first-aid kit and medical kit from the LH wardrobe, RH wardrobe, outboard galley, forward galley or aft galley, as applicable ([Figure 601](#)).

- (2) For aircraft equipped with 16 seats, remove the first-aid kit and medical kit from the LH forward stowage ([Figure 602](#)).

- (3) Do the inspection of the first-aid kit as follows:

- (a) Inspect the first-aid kit for violation and expiration date of the packaging of all the items of the kit.

NOTE: The items with the packaging violated must be replaced.

- (4) Do the inspection of the medical kit as follows:

- (a) Examine the medical kit to make sure that its permitted life is not expired and that the tell-tale wire is not broken.


NOTE: If the permitted life is expired or the tell-tale wire is broken, make sure that you can find all items in the kit and examine the date which shows the permitted life of each of them. Include or replace items, as applicable.

- (5) Put the first-aid kit and medical kit back on the LH wardrobe, RH wardrobe, outboard galley, forward galley or aft galley, as applicable ([Figure 601](#)).
- (6) For aircraft equipped with 16 seats, put the first-aid kit and medical kit back on the LH forward stowage ([Figure 602](#)).

EFFECTIVITY: ACFT MODEL(S) EMB-135

First-Aid Kit and Medical Kit - Inspection

Figure 601

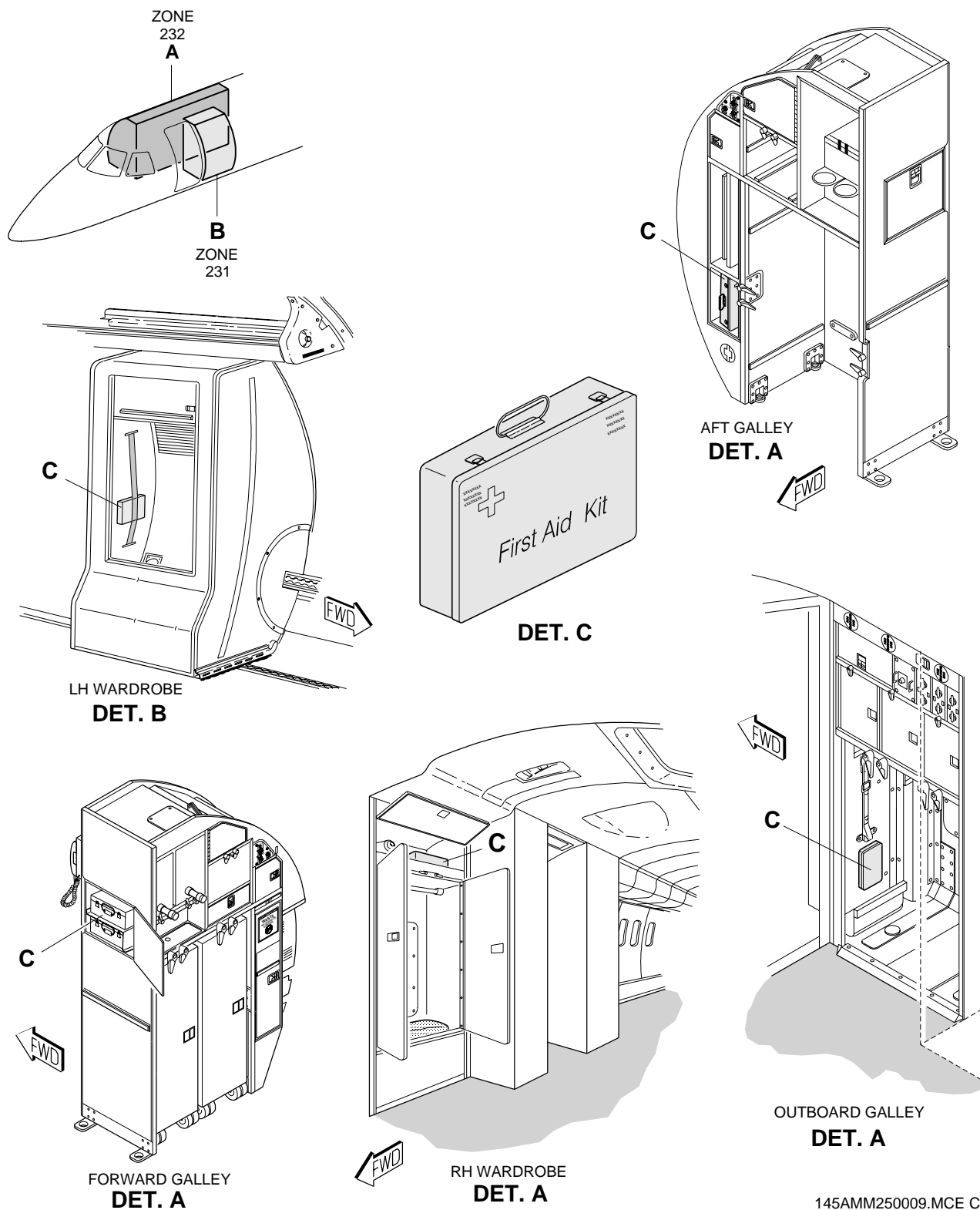


145AMM250009.MCE C

EFFECTIVITY: AIRCRAFT FOR 16 PASSENGERS

First-Aid Kit and Medical Kit - Inspection

Figure 602

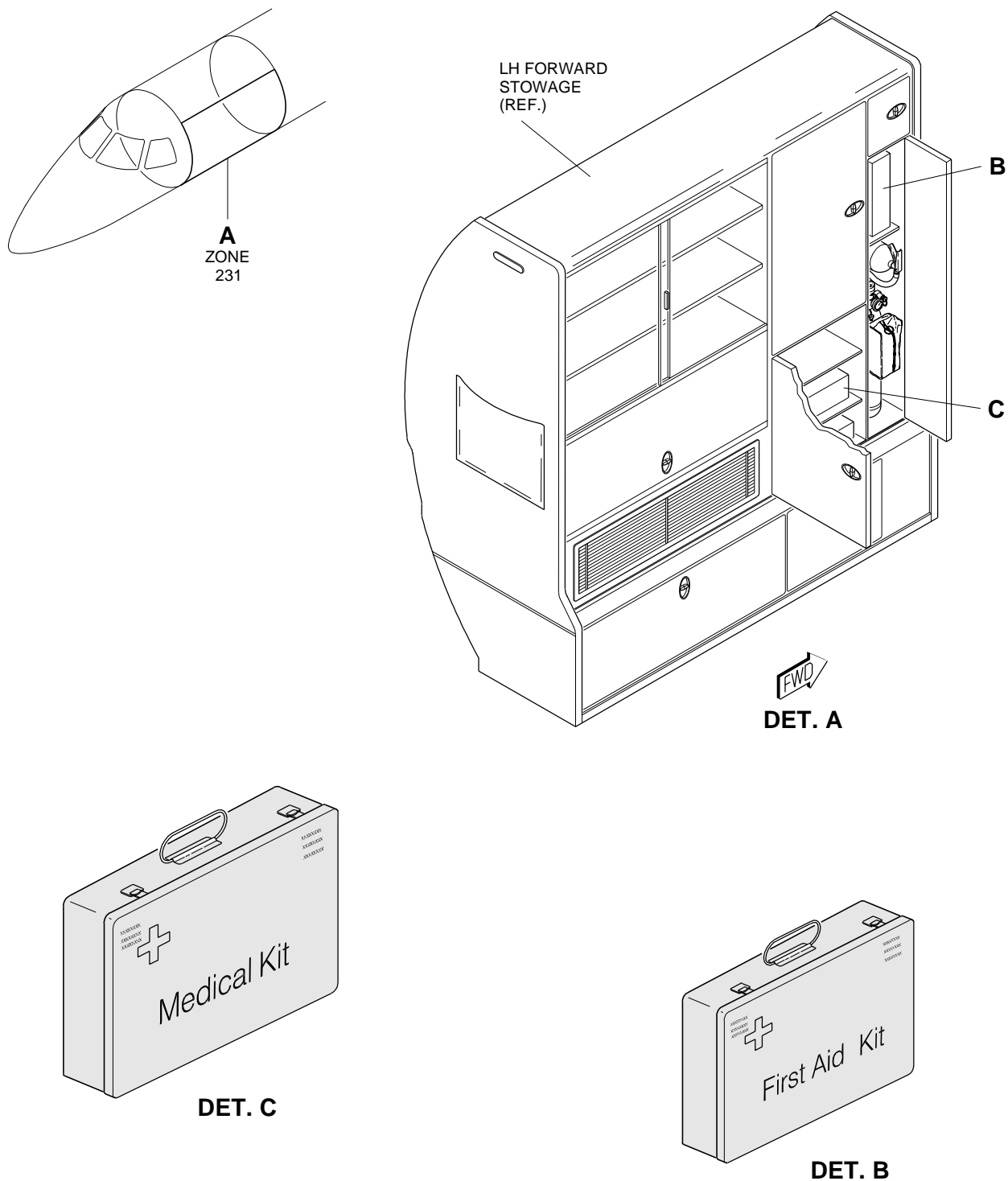

145AMM251656.MCE

